


# Leveraging Monitoring Technology at Point-of-Care

## BeneVision N-Series Monitors and Telemetry Pairing




Mindray optimizes essential monitoring at bedside with the pairing capability of the BeneVision N-Series and TM70 and TM80 telemeters, supporting patient mobility and migration from bedside to ambulation. By pairing the devices, patient transfer is seamless with no need to change electrodes, lead sets nor patient cable, creating an easy transition for both clinician and patient. When paired, the TM70/TM80 telemeter acts as an ECG, HR, SpO<sub>2</sub>, and respiration acquisition source, transmitting data to the N-Series bedside monitor, the BeneVision Distributed Monitoring System (DMS) and all deployed remote monitoring applications. This allows clinicians to easily flex patients from higher to lower acuity and back again, while maintaining a single patient sector on the BeneVision DMS, and the secure transmission of patient data to the EMR, ensuring a complete record.

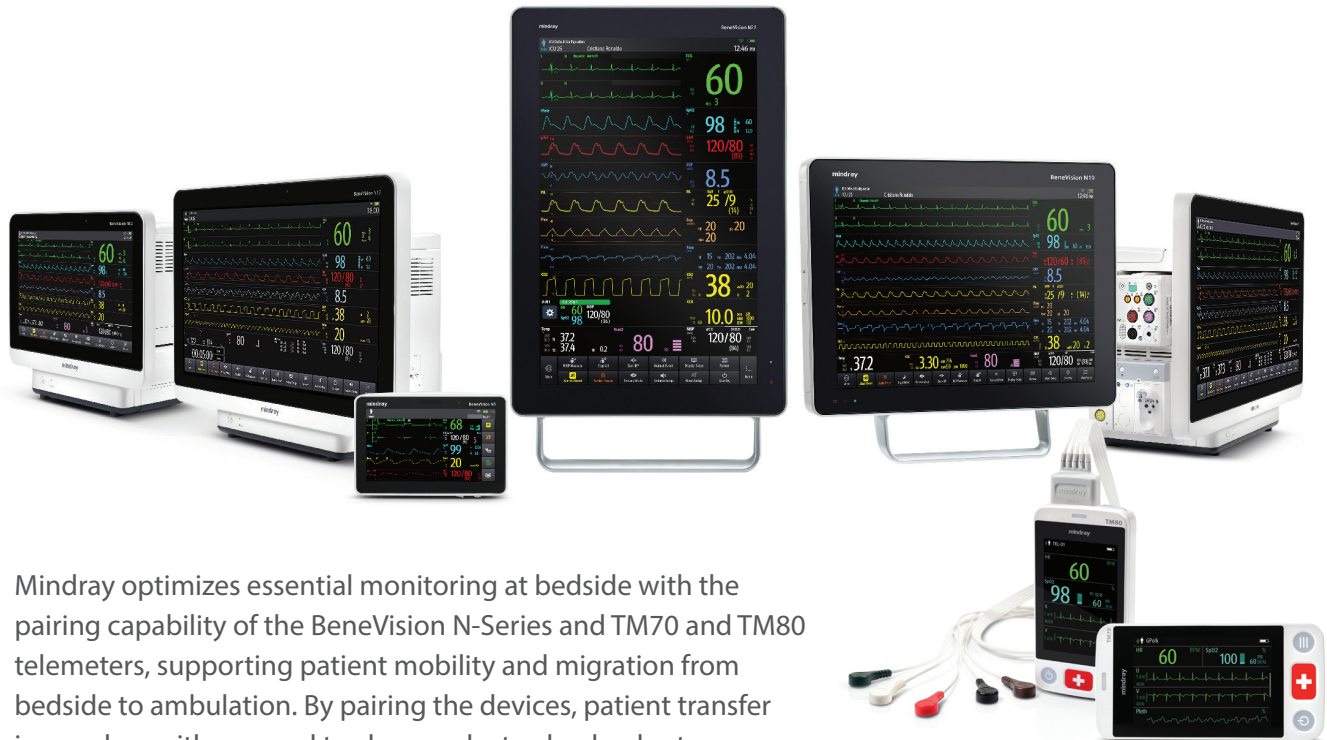

- Simplifies complex clinical workflow when transitioning patients in the ICU and stepdown
- Encourages patient mobility transitioning between rest and ambulation, supporting patient recovery
- Enables uninterrupted monitoring and a complete patient record for delivery of safe, quality care during device transition

# Leveraging Monitoring Technology at Point-of-Care

## BeneVision N-Series Monitors and Telemetry Pairing

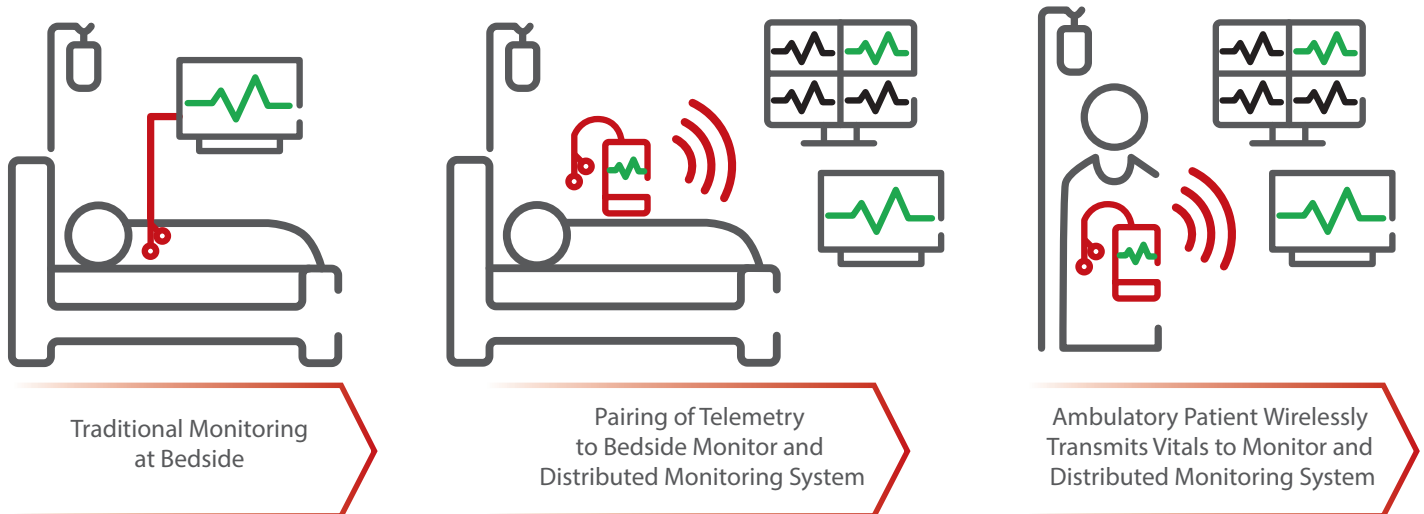

### Maximize Time Spent Caring for Your Patients

- With one touch on the N-Series monitor, access the Bedside Devices menu to look up telemetry devices for quick and efficient pairing.
- Initiate telemetry and monitor pairing at bedside, reducing disruption to patient care.
- Telemetry data is clearly displayed at the paired N-Series and on a shared patient sector at the BeneVision DMS providing flexible data access.
- Easily unpair the telemetry on the N-Series at bedside, saving valuable time in the process.

### Software Compatibility for Pairing

Mindray Device	Minimum Software Version Required
BeneVision N12, N15, N17, N19, N22	01.11.00
BeneVision TM70/TM80 with Alarms	4.01
BeneVision DMS	04.04.00
eGateway	7.3



Mindray North America  
 800 MacArthur Blvd., Mahwah, NJ 07430  
 Tel: 800.288.2121 Tel: 201.995.8000 Fax: 800.926.4275 [www.mindray.com](http://www.mindray.com)

©2021 Mindray DS USA, Inc.  
 Mindray® is a registered trademark of Shenzhen Mindray Bio-Medical Electronics Co., Ltd. All rights reserved.  
 All brands and product names are trademarks of their respective owners. 11/21 P/N: 0002-08-8163 Rev A

**mindray**  
 NORTH AMERICA