



ANESTHESIA WORKSTATION

Maximize Patient Safety

- Electronic vaporizers for safe, precise and reliable agent delivery
- High Flow Nasal Cannula (HFNC) helps clinicians intubate more easily
- Advanced ventilation modes to meet challenging patient demands (VCV, PCV, PCV-VG, CPAP/PS, SIMV modes standard with APRV optional)
- Volume Exchanger (VE) breathing system with visual indicator for ICU-level ventilation
- Protective ventilation toolkits to reduce incidence of post-operative pulmonary complications

Complement Clinical Workflow

- Intuitive user interface on rotatable 18.5" touchscreen monitor, with illuminated APL valve in manual mode
- Comprehensive, fast and simple system check, completed in 3.5 minutes, saves preparation time
- Save up to 10 customizable profiles
- Easy-to-read status dashboard enables quick and efficient troubleshooting

Enhance Efficiency

- Gas reduction tools (Optimizer, Flow Pause, e-AGSS, Agent Usage) provide both environmental and economic benefits
- Highly flexible integration options
- Safe and seamless connectivity and integration to meet diverse clinical needs
- Best-in-class 3-year warranty on A9 platform




The A9 anesthesia workstation has been designed around clinicians, with deep insights into clinical workflow, aiding in reduced clinical workload and maximized safety all-around. Mindray's innovative A9 platform retains traditional features while introducing electronic vaporizers to ensure the safety of patients throughout the perioperative period, complement the way clinicians work best and increase efficiency in the operating room.

mindray



ANESTHESIA WORKSTATION

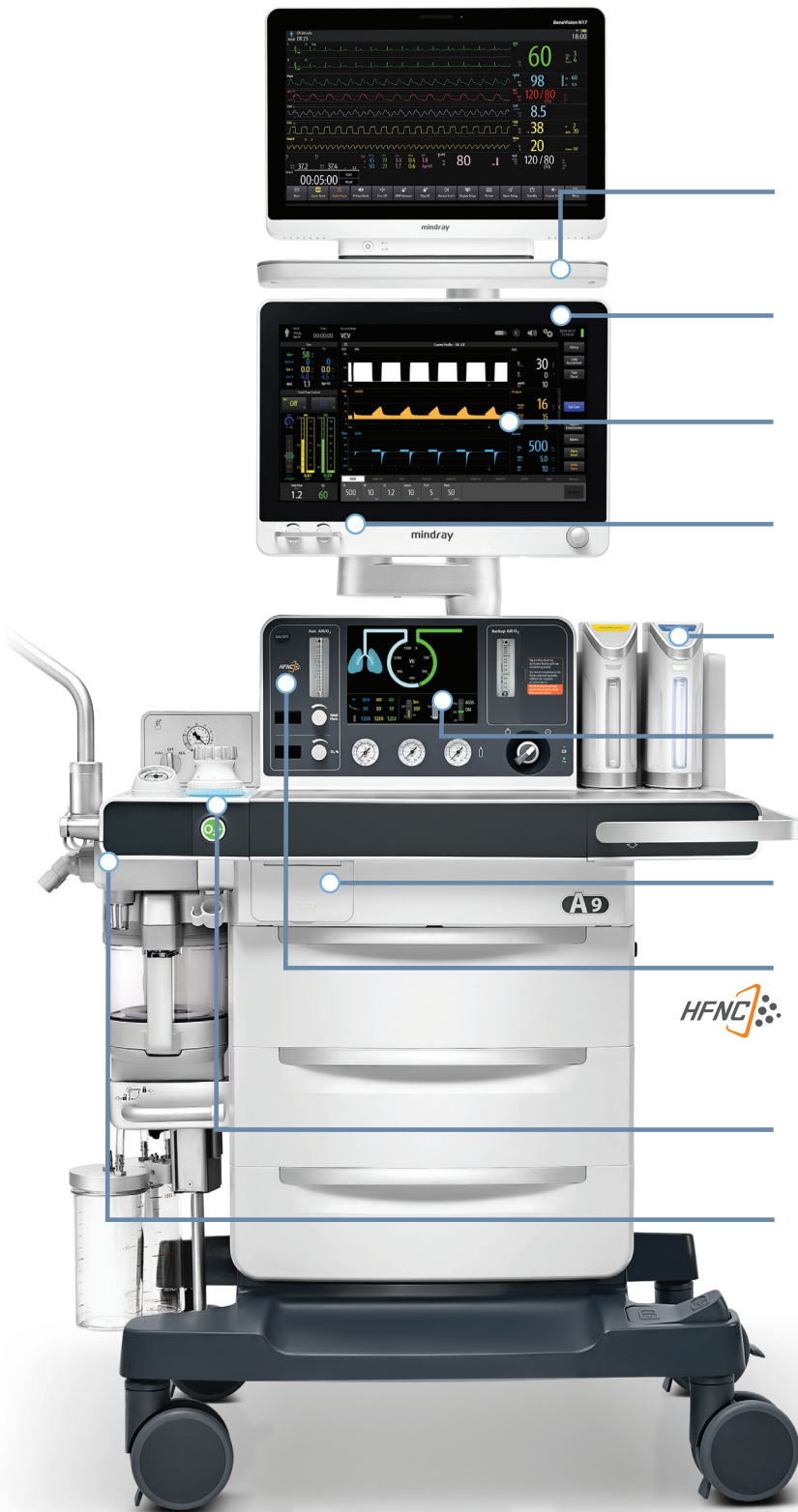

Highly flexible integration options work together with a variety of devices to save space in the operating room

18.5 inch capacitive touchscreen, with flexible rotation for 360 degrees of view

Powerful toolkits support confident decision-making for lung-protective ventilation

Electronic flowmeter with traditional easy-to-use knobs supports multiple setting methods

Precise electronic vaporizers, in a familiar position, notify user when agent gets low (A9 only)

VE indicator and status screen allow for quick visualization of the patient and system status to aid in rapid troubleshooting

ACGO, for use with non-rebreathing circuits, allows for changes in workflow



HFNC helps maintain safe oxygen saturation in patients to help clinicians intubate more easily, ideal for patients with poor oxygen saturation such as bariatric, pediatric, critically ill or those with a difficult airway.

Illumination around the APL valve in manual mode enables easy recognition of current operating mode

Integrated breathing circuit delivers precise ventilation



Mindray North America
800 MacArthur Blvd., Mahwah, NJ 07430
Tel: 800.288.2121 Tel: 201.995.8000 Fax: 800.926.4275 www.mindray.com

©2021 Mindray DS USA, Inc.
Mindray® is a registered trademark of Shenzhen Mindray Bio-Medical Electronics Co., Ltd. All rights reserved. All brands and product names are trademarks of their respective owners. Specifications subject to change without notice. Contact your Mindray Representative for the most current information. 07/21 P/N: 0002-08-1732 Rev A

