

BeneLink

CONNECTIVITY MODULE

Features and Benefits

- Simple, cost-effective solution to interface external medical devices to Mindray bedside patient monitors
- Compatible with Mindray BeneVision N-Series (N12, N15, N17, N19, N22) and Passport Series (12m, 17m) modular patient monitors
- Seamlessly integrates real-time data from a variety of ventilators, anesthesia systems and standalone devices
- BeneLink Devices Integration Window (on the Mindray bedside monitor, BeneVision DMS WorkStation or BeneVision CMS Viewer) displays integrated parameter data, settings and alarms
- Data displayed includes device type, up to 24 numerics, alarm status and device mode (i.e. ventilation)
- Parameters received from external medical devices are also included in the bedside monitor graphic and tabular trends
- Single-width module with four RJ-45 ports supports connectivity of up to four devices simultaneously
- All data collected through the BeneLink module securely passes to the EMR through the Mindray eGateway and BeneVision DMS (if applicable), providing a comprehensive patient record



Powerful Device Integration at the Bedside

The **BeneLink Connectivity Module** consolidates essential patient data at the primary patient monitor supporting a more complete patient assessment at the bedside. BeneLink is designed for use with high acuity Mindray patient monitors in the OR and ICU environments. Connecting as many as four external medical devices, BeneLink is a cost-effective solution for integration of ventilators, anesthesia systems and other standalone devices.

The BeneLink Devices Integration Window allows the display of essential patient data from third-party devices at the Mindray bedside monitor, BeneVision Distributed Monitoring System (DMS) workstation or via BeneVision CMS Viewer. At each, parameter data from the external devices is incorporated into the graphic and tabular trends, further aiding clinical assessment and decision-making. Such consolidation enables critical patient data and alarms to be displayed, saved, recorded, printed and calculated at the bedside monitor, thereby supporting highly efficient workflow and patient management.

mindray

BeneLink CONNECTIVITY MODULE



A-Series Advantage



BeneLink



N12


N17

N1

N22

N19

N15

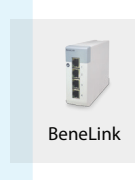

BeneLink



Ventilator

Supported Anesthesia Systems

Draeger	Apollo
	Fabius GS
	Fabius GS Premium (for BeneVision N-Series only)
	Fabius Tiro
	Perseus A500
GE	Primus (for Passport M-series monitors only)
	Aestiva 7900/7100
	Aespire 7900/7100
	AespireView (for BeneVision N-Series only)
	Aisys
Aisys	Aisys CS2 (for BeneVision N-Series only)
	Avance / Avance CS2
Maquet	FLOW-i
Mindray	A-Series (A7, A5, A4, A3)
	A8/A9 (for BeneVision N-Series only)

Supported Infusion Systems (BeneVision N-Series Monitors only)

B. Braun	Perfusor Space
Fresenius	Agilia Injectomat
	Agilia Injectomat MC
	Agilia Injectomat TIVA
	Agilia Volumat
	Agilia Volumat MC
	Fresenius Link+

Supported tcGas Monitoring Devices

Biometer	TCM40
	TCM 4 (for BeneVision N-Series only)
	TCM Tosca / TCM CombiM (for Passport M-Series only)
SenTec	SenTec Digital Monitor

Supported Ventilators

Carefusion	Avea (VOXP)(for BeneVision N-Series only)
	Vela
Draeger	Babylog 8000 / Babylog 8000 plus
	Evita 2 / Evita XL / Evita 2 Dura / Evita 4
	Infinity V500
	Babylog VN500
Savina	Savina 300
GE	Carescape R860
	Engstrom(Carestation)
Hamilton	C1/ C2/ C3 / T1 (Polling Protocol)(C1 and T1 are for BeneVision N-Series only)
	G5 (Block Protocol)
	G5 (Polling Protocol)
	Galileo (Polling Protocol)
Maquet	SERVO-i / SERVO-s
	SERVO-U / SERVO-N
Newport	E360
Philips	Respironics V60
Puritan Bennett	PB840 (SNDF protocol)
	PB840 (SNDA protocol)
	PB980 (SNDF protocol)
	PB980 (SNDA protocol)
ResMed	VSIII (for Passport M-Series only)

Supported Other Devices

Organon (NMT monitoring)	TOF Watch SX*
---------------------------	---------------

Mindray North America

800 MacArthur Blvd., Mahwah, NJ USA 07430

Tel: 800.288.2121 or 201.995.8000 Fax: 800.926.4275 www.mindray.com

©2022 Mindray DS USA, Inc. Mindray® is a registered trademark of Shenzhen Mindray Bio-Medical Electronics Co., Ltd. All rights reserved. All brands and product names are trademarks of their respective owners. Specifications subject to change without notice. Contact your Mindray Representative for the most current information. 01/22 P/N: 0002-08-8130 Rev C

mindray
NORTH AMERICA